

Nordent®  
**InstRenew<sup>2</sup>**  
SHARPENING ASSISTANT

# Just Like New... In Seconds!



It's Like Having Another  
Assistant in the Office!  
Get Professional  
Results Everytime.

The Most Precise  
Sharpening  
Assistant is  
Faster &  
Easier...

▶ InstRenew<sup>2</sup>  
Does ALL  
the Work

▶ Gives You Exact Angles  
EVERYTIME

## THE NORDENT GUARANTEE

We guarantee, you will agree Nordent  
Instruments offer the highest quality  
and best performance...  
or your money back.


SSM2000  
(MSRP \$369.50)

**Nordent**  
MANUFACTURING INC.

[www.nordent.com](http://www.nordent.com)

06.06.15M.MK301

# Nordent® **InstRenew<sup>2</sup>** SHARPENING ASSISTANT

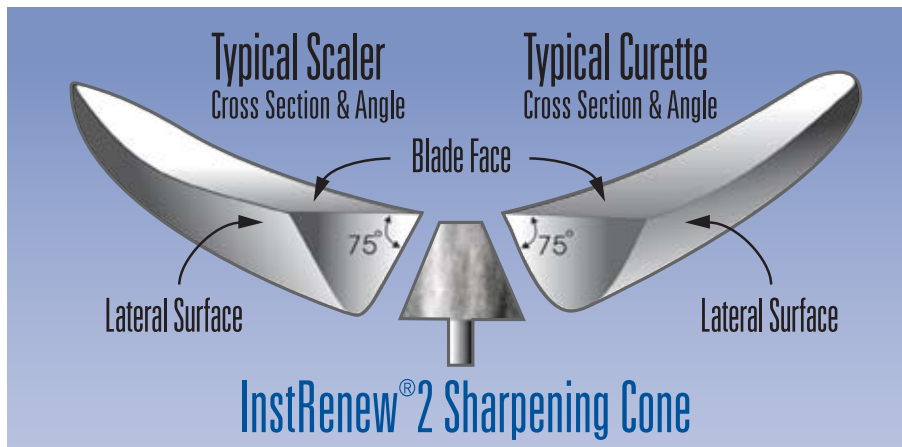

Every Brand and Every Pattern of scalers and curettes have one thing in common. The angle between the blade face and the lateral surface is always 75°.

## 3-Step Method

For Every Hygiene Instrument



### POSITION THE BLADE,

The Blade Positioner automatically locks **EVERY instrument** into the correct position for sharpening by holding the blade face horizontal.



### GRASP THE TIP,

The Tip Clamp allows you to easily grasp the instrument on any straight portion of the shank and maintain the correct position while you sharpen.

**That's It, You're Done!**



### SHARPEN THE BLADE!

The Sharpening Cone is cut at the correct 75° angle for **EVERY brand** and **EVERY pattern** of scalers and curettes. Simply trace each cutting edge over the sharpening cone.

**Nordent**  
MANUFACTURING INC.  
610 Bonnie Lane • Elk Grove Village, IL 60007  
Phone: 847.437.4780 • Toll Free: 800.966.7336  
[www.nordent.com](http://www.nordent.com)